

E4 XLS+ Multichannels

Precision Across All Channels



Easy operation

The E4 XLS+ works the way you think. Its color graphics, large screen and precision joystick have revolutionized electronic pipetting. The simple interface makes adjusting settings quick and easy.



Channel consistency

The LTS™ cylindrical tip design assures fast tip mounting, consistent sample loading and leak-free seals across all channels. In addition to LTS, the new push bar design in all XLS+ multichannel pipettes assures exceptionally even and consistent aspiration across all channels.



Simplify repetitive pipetting

Multi-dispense takes the drudgery out of repetitive pipetting by calculating how many aliquots can be dispensed from a single pickup. Autopace simplifies repeat dispensing – click the joystick once and repeat dispense across a plate.



Save your protocols

Store multiple custom protocols in the E4's on-board memory. Simultaneously review and select your protocols with a simple movement of the joystick. No more re-entering protocols over and over again!



E4 XLS+ Electronic Multichannels

Speed Your Workflow


The E4™ XLS+™ multichannel pipette is a perfect tool for high-throughput work. Fast loading, lightweight, easy to handle and precision and consistency across all channels. Pipette settings, protocols and service alarms can be password protected for GLP/GMP compliance. What's more, GLP data, such as service records, cycles and status data is completely tamper proof.

E4 XLS+ Electronic Multichannel Pipette

Precise, High-throughput Pipetting



- LTS™ for consistent sample loading across all channels and up to 80% lower tip ejection forces
- Standard modes include Basic pipetting, Multi-dispense and Mix
- Advanced modes include sequential volume programming for complex protocols, titration, dilution and reverse mode pipetting
- Control of aspirate, dispense and mix speeds for challenging liquids
- Optional blowout and automatic dispense settings for precision pipetting
- Password protect pipette settings for GLP security
- Save your favorite protocols
- Simple or complex: configure the E4 for any task
- FPLC in your tip: on board PureSpeed™ sample preparation protocol
- Multilingual display and on-screen Help
- Easily accessible service records and alarms
- RFID enabled for advanced asset management
- Optional Rapid Charge Stand charges three pipettes in one hour

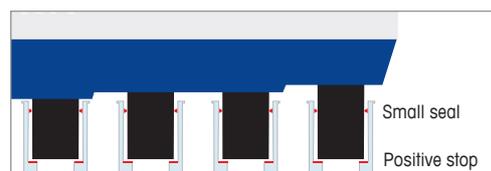

LTS™ LiteTouch™ Tip Ejection System

- Small seal and positive stop means consistent sample loading – across all channels, every time!
- No nozzle o-rings
- Lower forces – less fatigue
- Progressive tip ejector

Technical Data

Specifications

| Model | Volume (µL) | Increment (µL) | Accuracy | | Precision | |
|---------|-------------|----------------|----------|--------|-----------|--------|
| | | | % | µL (±) | % | µL (≤) |
| 10 µL | 1,0 | 0,01 | 4 | 0,04 | 3 | 0,03 |
| | 5,0 | | 1,5 | 0,075 | 0,6 | 0,03 |
| | 10,0 | | 1,0 | 0,1 | 0,4 | 0,04 |
| 20 µL | 2 | 0,02 | 7,5 | 0,15 | 2,0 | 0,04 |
| | 10 | | 1,5 | 0,15 | 0,5 | 0,05 |
| | 20 | | 1,0 | 0,2 | 0,3 | 0,06 |
| 50 µL | 5 | 0,05 | 3,5 | 0,18 | 1,5 | 0,075 |
| | 25 | | 1,2 | 0,3 | 0,4 | 0,1 |
| | 50 | | 0,8 | 0,4 | 0,2 | 0,1 |
| 200 µL | 20 | 0,2 | 2,5 | 0,5 | 1,0 | 0,2 |
| | 100 | | 0,8 | 0,8 | 0,25 | 0,25 |
| | 200 | | 0,8 | 1,6 | 0,15 | 0,3 |
| 300 µL | 30 | 0,2 | 2,5 | 0,75 | 1,0 | 0,3 |
| | 150 | | 0,8 | 1,2 | 0,25 | 0,375 |
| | 300 | | 0,8 | 2,4 | 0,15 | 0,45 |
| 1200 µL | 100 | 1 | 3,6 | 3,6 | 0,6 | 0,6 |
| | 600 | | 0,8 | 4,8 | 0,2 | 1,2 |
| | 1200 | | 0,8 | 9,6 | 0,15 | 1,8 |

Rainin, LTS, LiteTouch, E4, XLS, XLS+ and Pipetting 360° are trademarks of Rainin Instrument, LLC. LabX and Direct Pipette-Scan are trademarks of METTLER TOLEDO.

Ordering Information

| Cat. No. | MT-Order No. | Volume Range |
|-------------------|--------------|--------------|
| 8-Channel | | |
| E8-10XLS+ | 17013792 | 0,5–10 µL |
| E8-20XLS+ | 17013793 | 2–20 µL |
| E8-50XLS+ | 17013794 | 5–50 µL |
| E8-200XLS+ | 17013795 | 20–200 µL |
| E8-300XLS+ | 17013796 | 20–300 µL |
| E8-1200XLS+ | 17014498 | 100–1200 µL |
| 12-Channel | | |
| E12-10XLS+ | 17013797 | 0,5–10 µL |
| E12-20XLS+ | 17013798 | 2–20 µL |
| E12-50XLS+ | 17013799 | 5–50 µL |
| E12-200XLS+ | 17013800 | 20–200 µL |
| E12-300XLS+ | 17013801 | 20–300 µL |
| E12-1200XLS+ | 17014499 | 100–1200 µL |

Accessories

| Cat. No. | MT-Order No. | Description |
|------------|--------------|---|
| E4-RCS | 17012332 | Rapid Charge Stand |
| RFID-KIT | 17011966 | RFID Kit including Rainin RFID Reader and LabX Direct Pipette-Scan Software |
| RFID-RDR | 17011964 | Rainin RFID Reader |
| LABX-PIPET | 17011965 | LabX Direct Pipette-Scan Software |

Note: Li-Ion battery and Wall Power Supply included.



Rainin Instrument, LLC

7500 Edgewater Drive, Oakland, CA 94621
Phone +1 510 564 1600

a METTLER TOLEDO Company

Subject to technical changes
© 10/2014 Rainin Instrument, LLC
Printed on demand from online content. H17700764 (Rev G)

www.mt.com/rainin

For more information